



Shoulder System

UNITED Oncology System.

A COMPLETE SOLUTION TO YOUR PATIENTS



Hip System



Knee System(Distal Femur)



Knee System(Proximal Tibia)

Preface

The UNITED Oncology joint system has been developed for more than 15 years. Through the modification done by the manufacturers and major designers, the major design of these prostheses is completely competitive to others in the market. The purpose of designing these prostheses is to meet the requirement of those patients who need a special prosthesis for their particular indication. Therefore, being able to provide fast delivery as well as custom-made shape and size for individual patient is most important. The UNITED Oncology System is based on modular combinations of the implants that keep most of components on shaft, and if, it is necessary to customize the critical parts to meet different patients' situation.

The first patient underwent surgery in 1996. Until now, there have been over seven hundred patients who benefited from the UNITED Oncology System for the actual sizing and efficient delivery.


With fast custom ability of the UNITED Oncology System, surgeons can adapt in their ideas or modifications for the patients who need special attention in a very short delivery time. The modular feature of the prostheses allows surgeons to choose various components depending upon the need of patient and full interchangeability, thus makes no surprise during the surgery.

In order to reduce the osteolysis, the polyethylene of the tibia insert and locking bushing were both increased in thickness to enhance resistance to polyethylene wear. The modularity, maximum range of motion, patented proprietary locking mechanism, are some of the advantages that the UNITED Oncology System provides to make the surgery much more complied to the surgeons' requirements.

Precise instruments to conduct an accuracy surgery are of much important to the UNITED Oncology System. The correct design of instruments can help surgeons to make the surgical procedure more accurate, simpler and faster. To obtain the minimum bone cut, to keep long term stable fixation, and special instruments for special requirement are all considerable in the design of the system.




UNITED Oncology Hip System

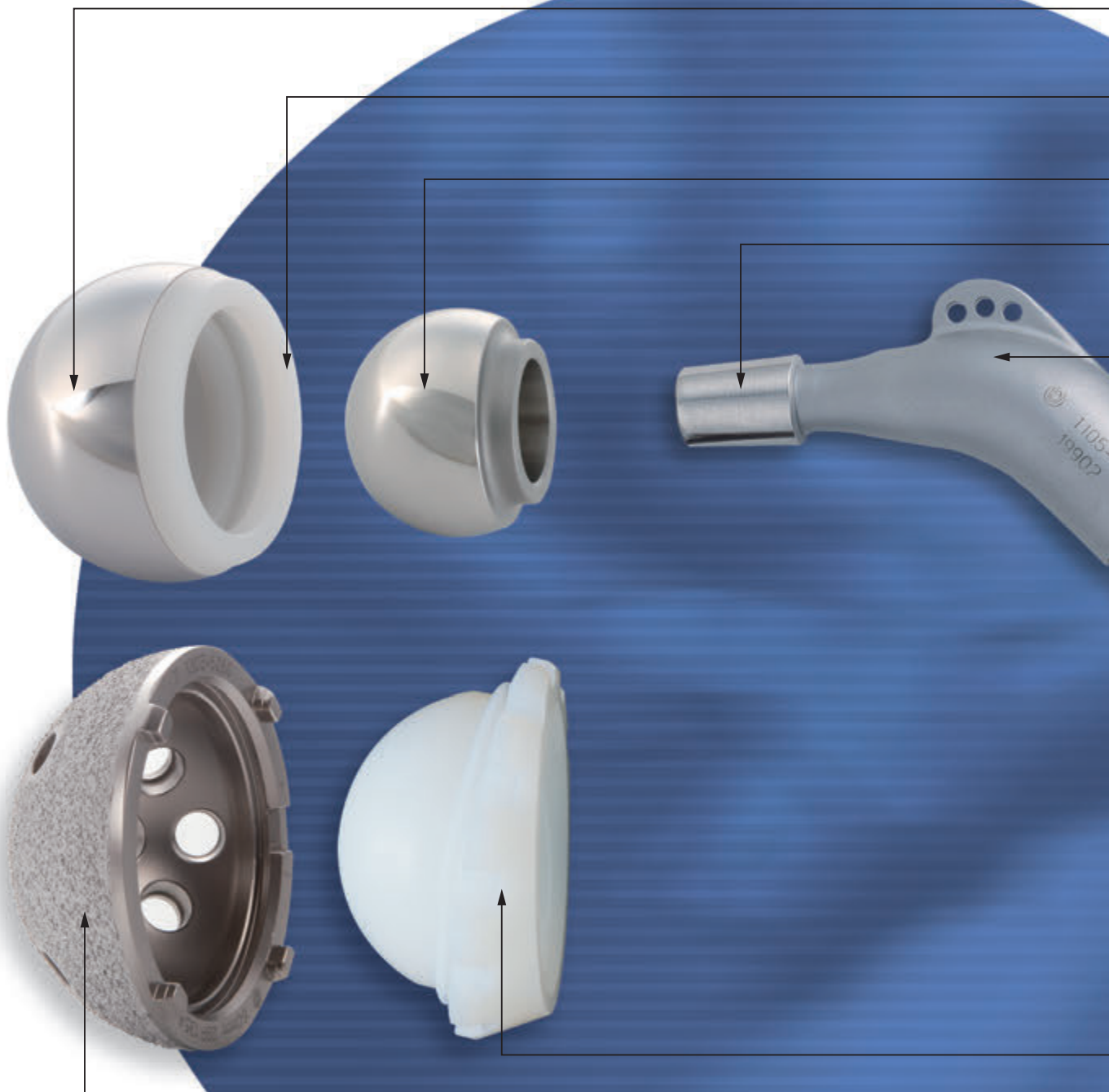

Case 1:
male, 40 yr, right proximal femur severe
bone deficiency secondary AVN, post-op 3 yr.



Case 2:
female, 65 yr, ovarian cancer, right proximal
femur metastasis, post-op 1 yr.

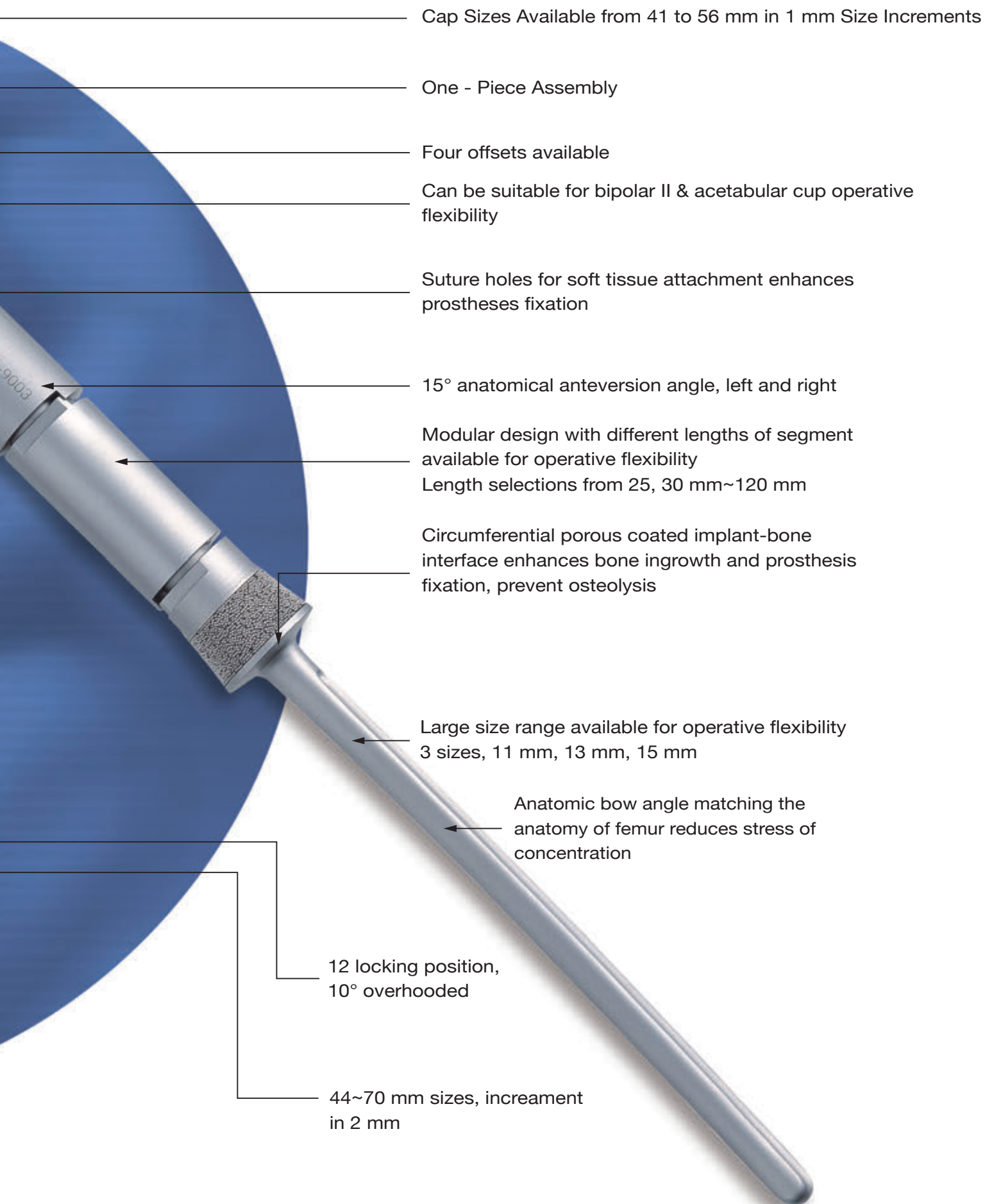

UNITED Oncology Knee System (Distal Femur)

Large size range available for operative flexibility

Anatomic bow angle matching the anatomy of femur
reduces stress concentration
3 sizes, 11 mm, 13 mm, 15 mm

Circumferential porous coated implant-bone interface
enhances bone ingrowth and prosthesis fixation,
prevents osteolysis


Modular design with different lengths of segment available
for operative flexibility
Length selections from 25, 30 mm~120 mm


Raised anterior slope prevents hyperextension

Locking mechanism design makes assembly easier

Spherical contact surface reduces stress concentration,
specific stop mechanism allows $\pm 20^\circ$ rotational freedom

Three different sizes selection, for maximum bone coverage

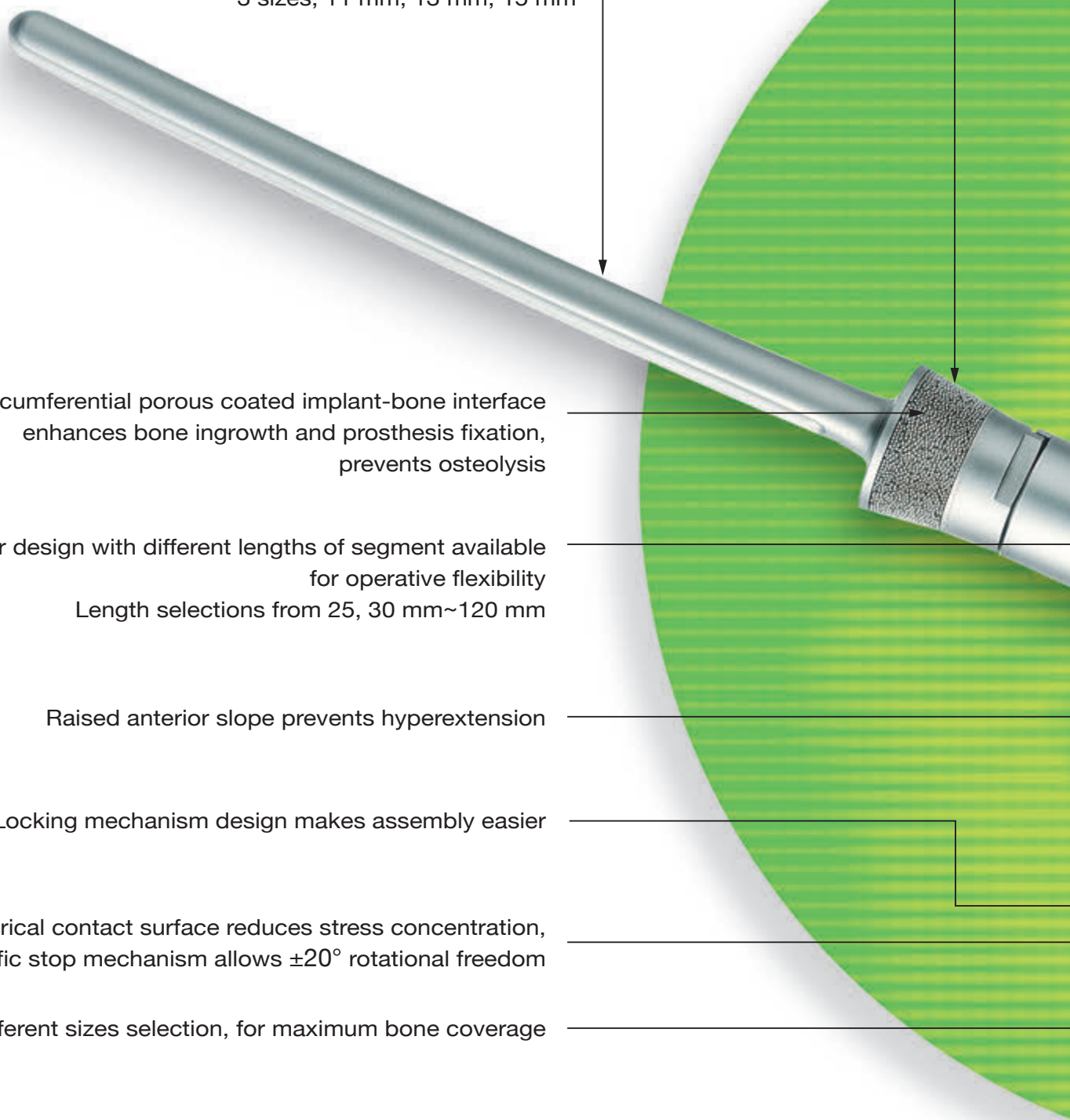

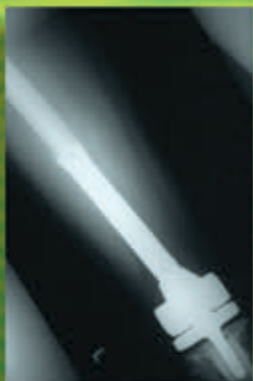

Case 1:
male, 17 yr, malignant
fibrosarcoma, post-op 3 yr.




Case 2:
male, 13 yr, malignant osteosarcoma, post-op 3yr.



Case 3:
male, 44 yr, MFH, left-distal femur, post-op 1yr.




UNITED Oncology Knee System (Proximal Tibial)

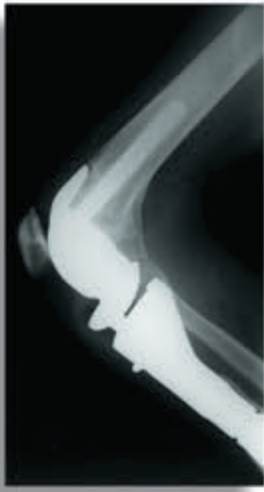


Case 1:
female, 16 yr, malignant
Grant-celltumor, post-op 3 yr.



Case 2:
male, 14yr, malignant
fibrosacoma, post-op 2 yr.

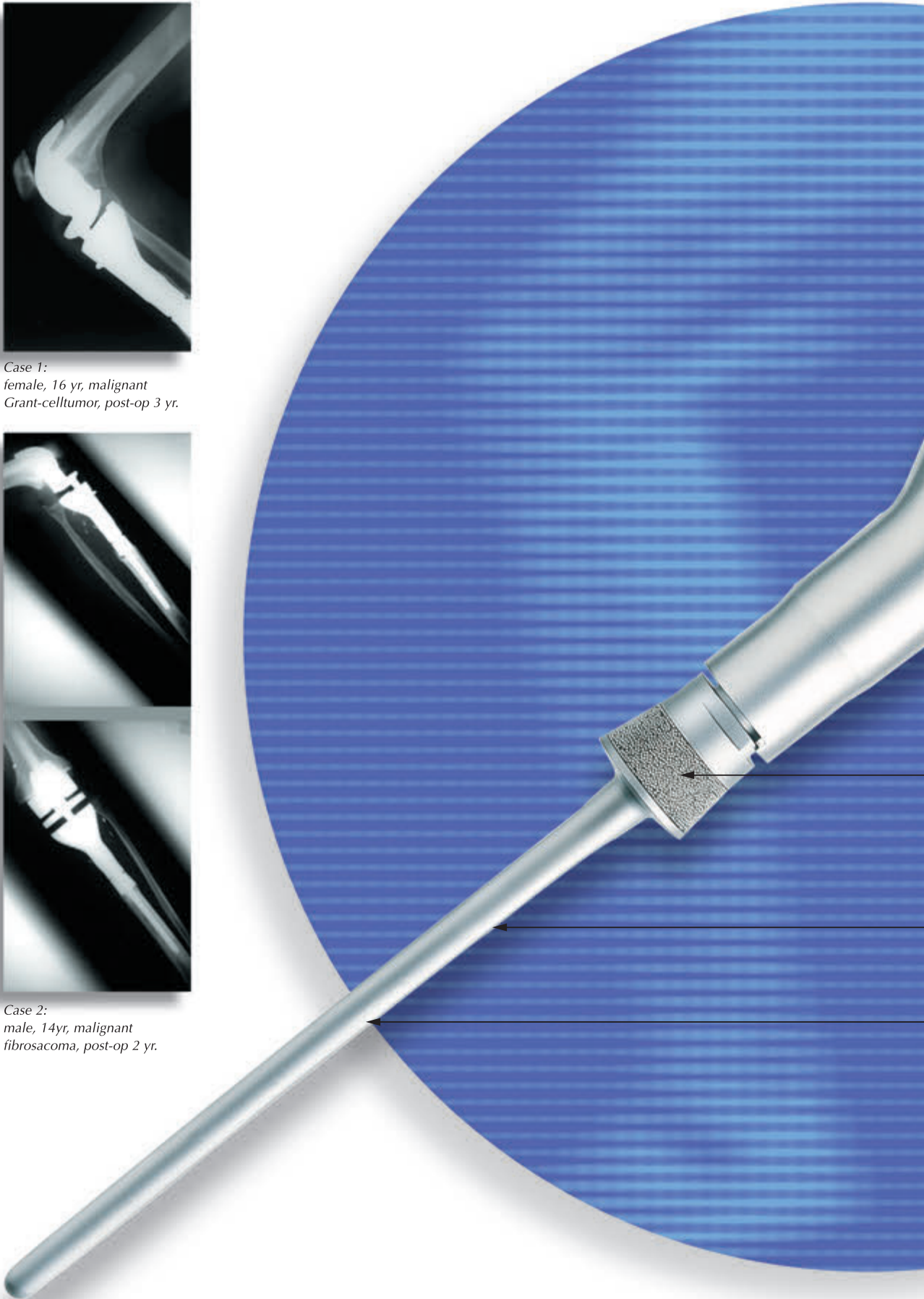

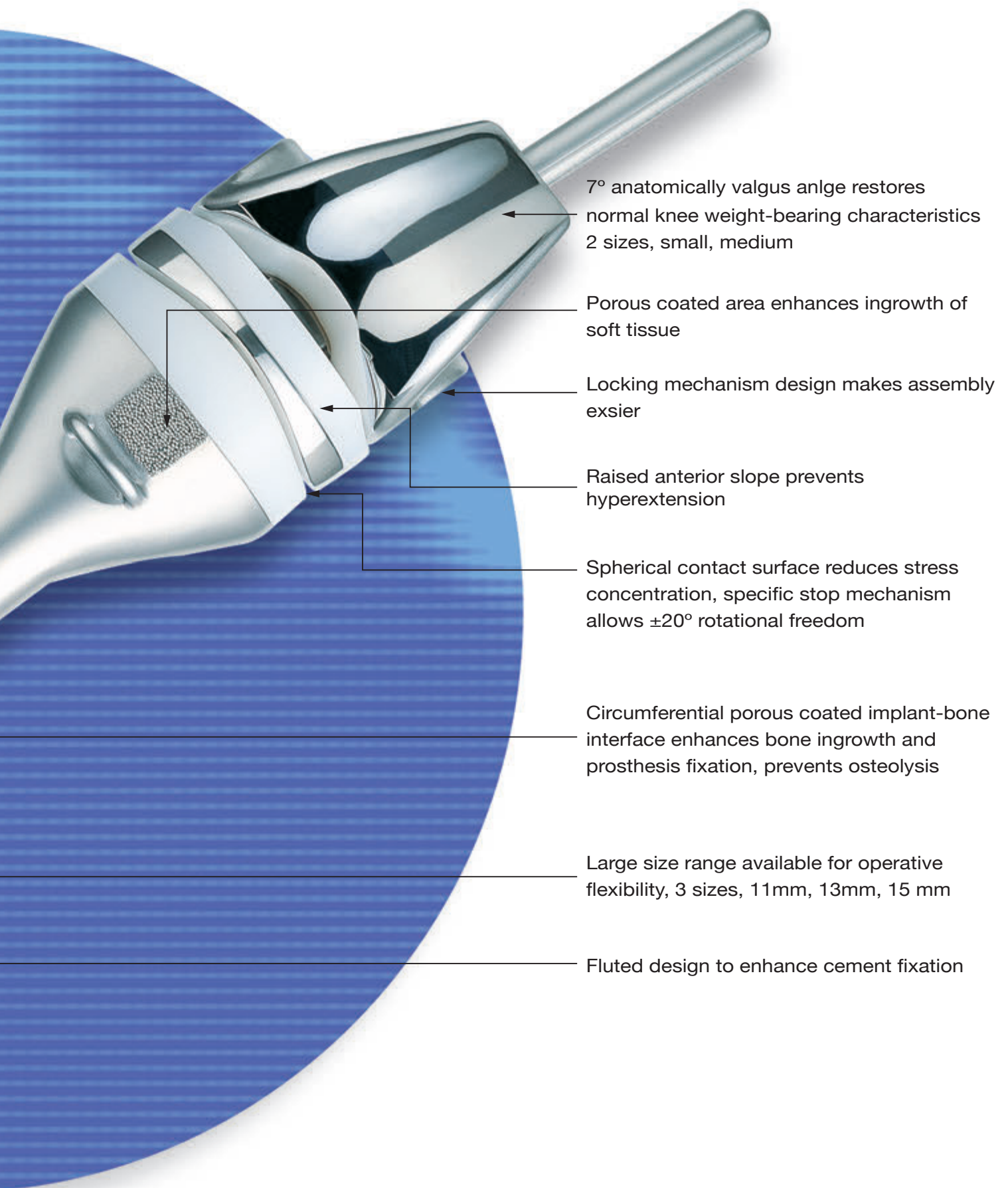

7° anatomically valgus angle restores normal knee weight-bearing characteristics
2 sizes, small, medium

Porous coated area enhances ingrowth of soft tissue

Locking mechanism design makes assembly easier

Raised anterior slope prevents hyperextension

Spherical contact surface reduces stress concentration, specific stop mechanism allows $\pm 20^\circ$ rotational freedom

Circumferential porous coated implant-bone interface enhances bone ingrowth and prosthesis fixation, prevents osteolysis

Large size range available for operative flexibility, 3 sizes, 11mm, 13mm, 15 mm

Fluted design to enhance cement fixation

UNITED Oncology Shoulder System



Three different sizes, 40 mm, 42 mm, 44 mm for optimum fixation of the patients.

Suture holes for soft tissue attachment enhances prostheses fixation

Modular design with different lengths of segment available for operative flexibility
Length selections from 25~80 mm

Circumferential porous coated implant-bone interface enhances bone ingrowth and prosthesis fixation, prevents osteolysis

Large size range available for operative flexibility
2 sizes, 7 mm, 9 mm


Fluted design to enhance cement fixation




Case 1:
female, 18 yr, malignant
condrosarcoma, post-op 1 yr.



Case 2:
male, 15 yr, MFH, post-op 3 yr.

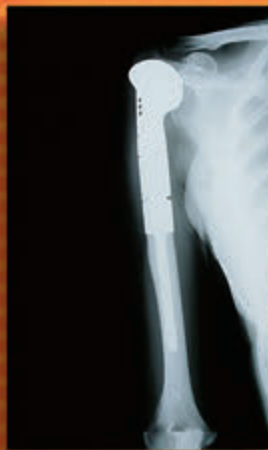

Case 3:
female, 16 yr, malignant
osteosarcoma, post-op 2 yr.

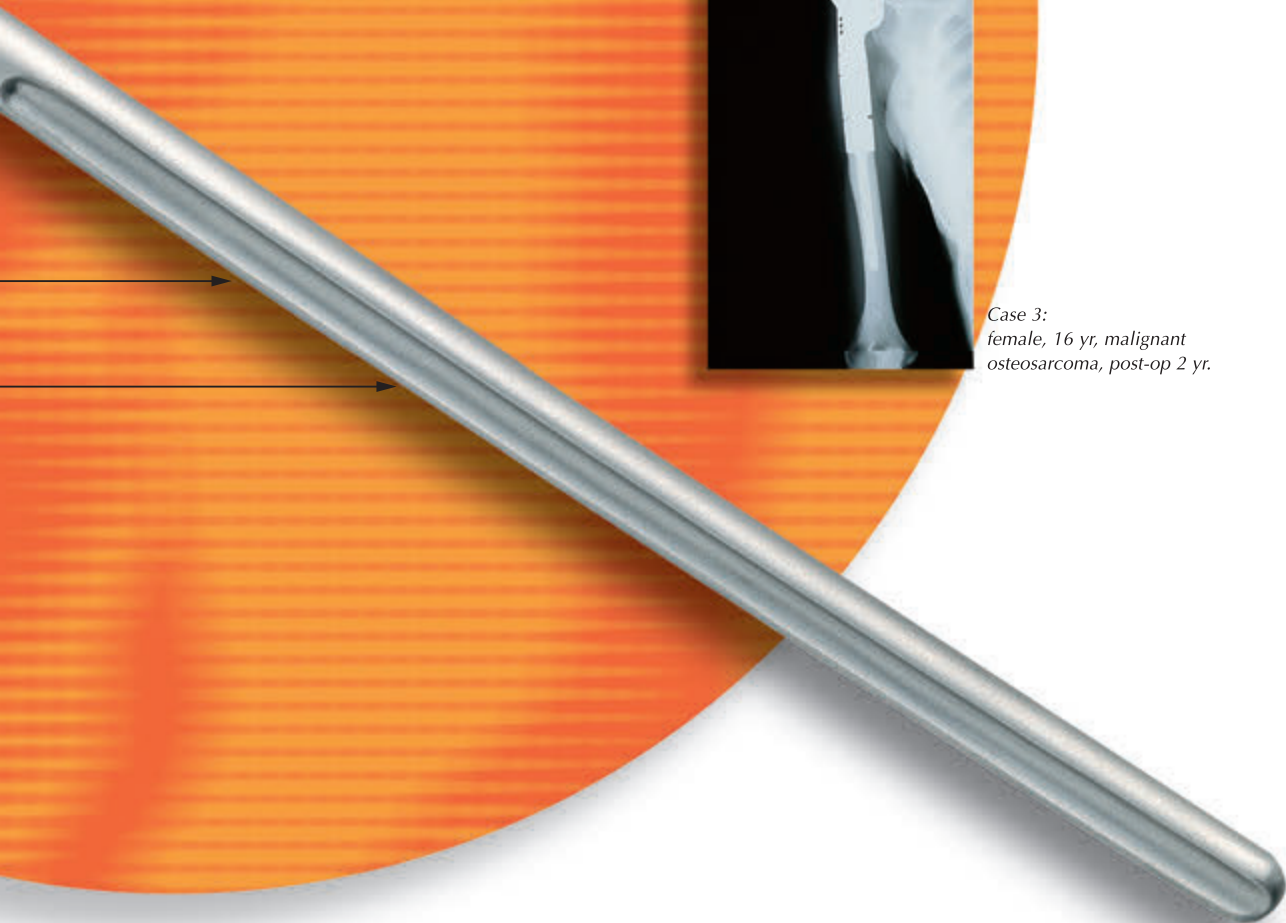



Each Step We Care

Contact Us

■ Taiwan

United Orthopedic Corporation

12F, No.80, Sec.1, Chenggong Rd., Yonghe Dist.,
New Taipei City 23452, Taiwan
Tel: +886 2 2929 4567 Fax: +886 2 2922 4567

■ China

United Medical Instrument Co., Ltd.

Room 402, Building 7, No.697, Ling Shi Road,
ShangHai 200072, China
Tel: +86 21 3661 8055 Fax: +86 21 6628 3031

■ Japan

United Biomech Japan Inc.

Ginyo Bldg 4th floor, 2-9-40 Kitasaiwai, Nishi-Ku
Yokohama City, Kanagawa, 220-0004, Japan
Tel: +81 45 620 0741 Fax: +81 45 620 0742

■ USA

UOC USA Inc.

20, Fairbanks, Suite 173, Irvine 92618, CA
Tel: +1 949 328 3366 Fax: +1 949 328 3367

■ EU Representative

mdi Europa GmbH

Langenhagener Strasse 71, Langenhagen 30855, Germany
Tel: +49 511 3908 9530 Fax: +49 511 3908 9539

■ Europe

UOC (Suisse) SA

Avenue Général Guisan 60A, Pully 1009, Switzerland
Tel: +41 21 634 70 70 Fax: +41 21 634 70 76

UOC (France)

21, Rue de la Ravinelle, Nancy 54000, France
Tel: +33 3 83 23 39 72 Fax: +33 3 83 23 39 10



www.uoc.com.tw