



RepliCase[™] • Hip Cement Spacer Mold

An Instrument Set to Reproduce Hip Cement Spacer with Ease

- Total hip type spacer prosthesis prevent dislocation
- 3 sizes are available
- Made by medical grade silicone to allow repeated autoclave



Surgical Procedures

1. Size Selection

Measured Acetabular Size	50, 52	54	56 or above
Suggested Size	S	M	L

2.




3.



4.



5.

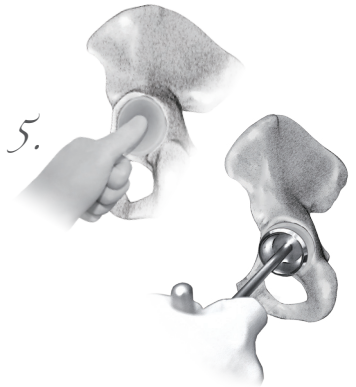


*Note: Bone cement and antibiotics are NOT provided.

For acetabular size under 50 mm, acetabular spacer cup might not be necessary, size selection should base on the head diameter. At least two batches of bone cement is needed for femoral spacer prosthesis.

The amount of bone cement required for spacer cup depends on acetabular size and bone defect.


Specification

Size	S	M	L
Cup Diameter (mm)	46	49	52
Head Diameter (mm)	42	45	48
Offset (mm)	36	39	43
Head Neck Angle	130°	130°	130°
Stem Length (mm)	145	145	145
Distal Diameter (mm)	9	10	11

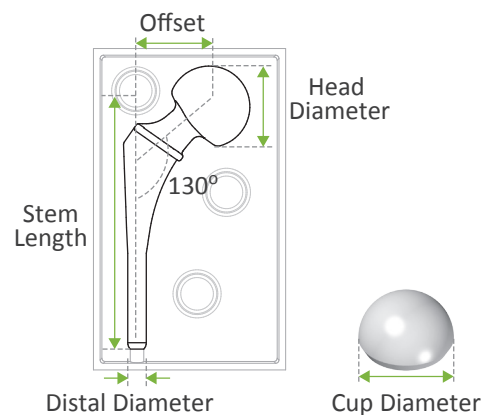

Order Information

RepliCase™ · Hip
Cement Spacer Mold

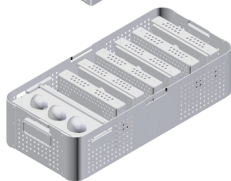

Catalog No. Description

9909-5009	S, 42 mm
9909-5010	M, 45 mm
9909-5011	L, 48 mm

Catalog No. Description

9909-8110	RepliCase Cement Mold Case
-----------	----------------------------

Catalog No. Description

9909-1101	Acetabular Cup Molder Handle
-----------	------------------------------

Catalog No. Description

9909-5046	Acetabular Cup Molder, 46 mm
9909-5049	Acetabular Cup Molder, 49 mm
9909-5052	Acetabular Cup Molder, 52 mm



Office

■ Taiwan

12F, No. 80, Sec. 1, Chenggong Rd, Yonghe Dist.,
New Taipei City 23452, Taiwan
Tel: +886 2 2929-4567 Fax: +886 2 2922-4567

■ USA

UOC USA Inc.
20 Fairbanks, Suite 173
Irvine, CA 92618
Tel: +1 949-328-3366 Fax: +1 949-328-3367

■ China

United Medical Instrument Co. Ltd.
Room 3-1-201, No.45, Jiangchang First Road, Shanghai, China
Tel: +86 21 3661-8055 Fax: +86 21 6628-3031

■ EU Representative

mdi Europa GmbH
Langenhagener Strasse 71, 30855 Langenhagen, Germany
Tel: +49-511-3908 9530 Fax: +49-511-3908 9539



www.uoc.com.tw

All Rights Reserved, ©2015 United Orthopedic Corp.

01CS-B-01-01