



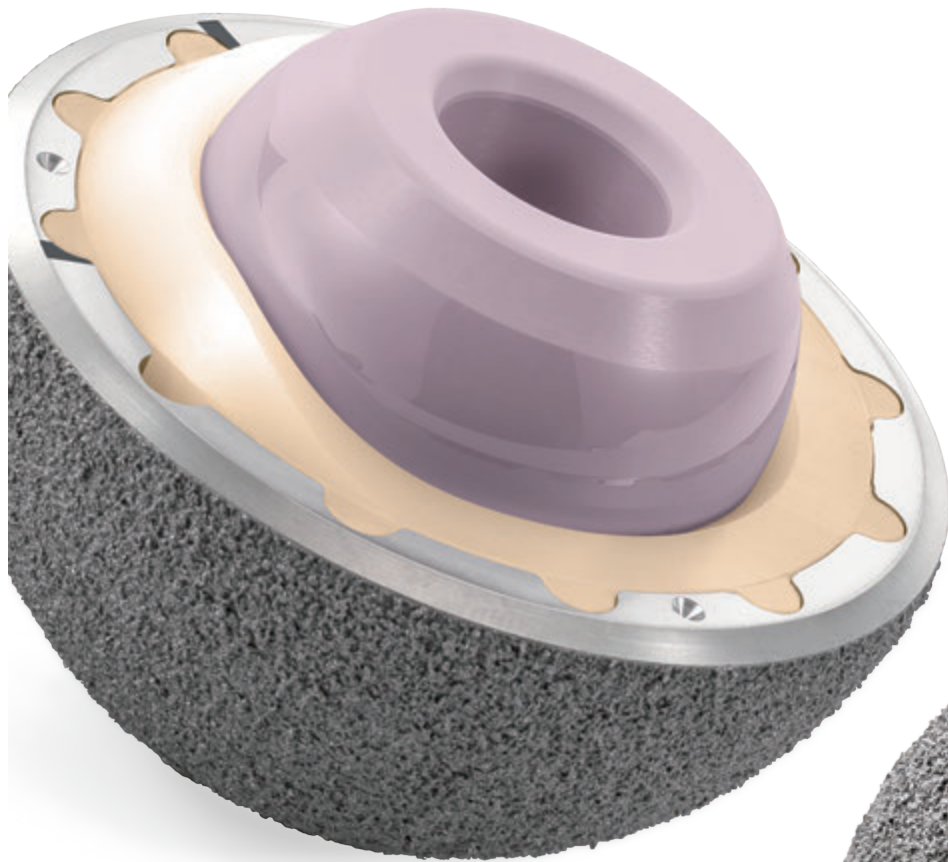


U-Motion II⁺™

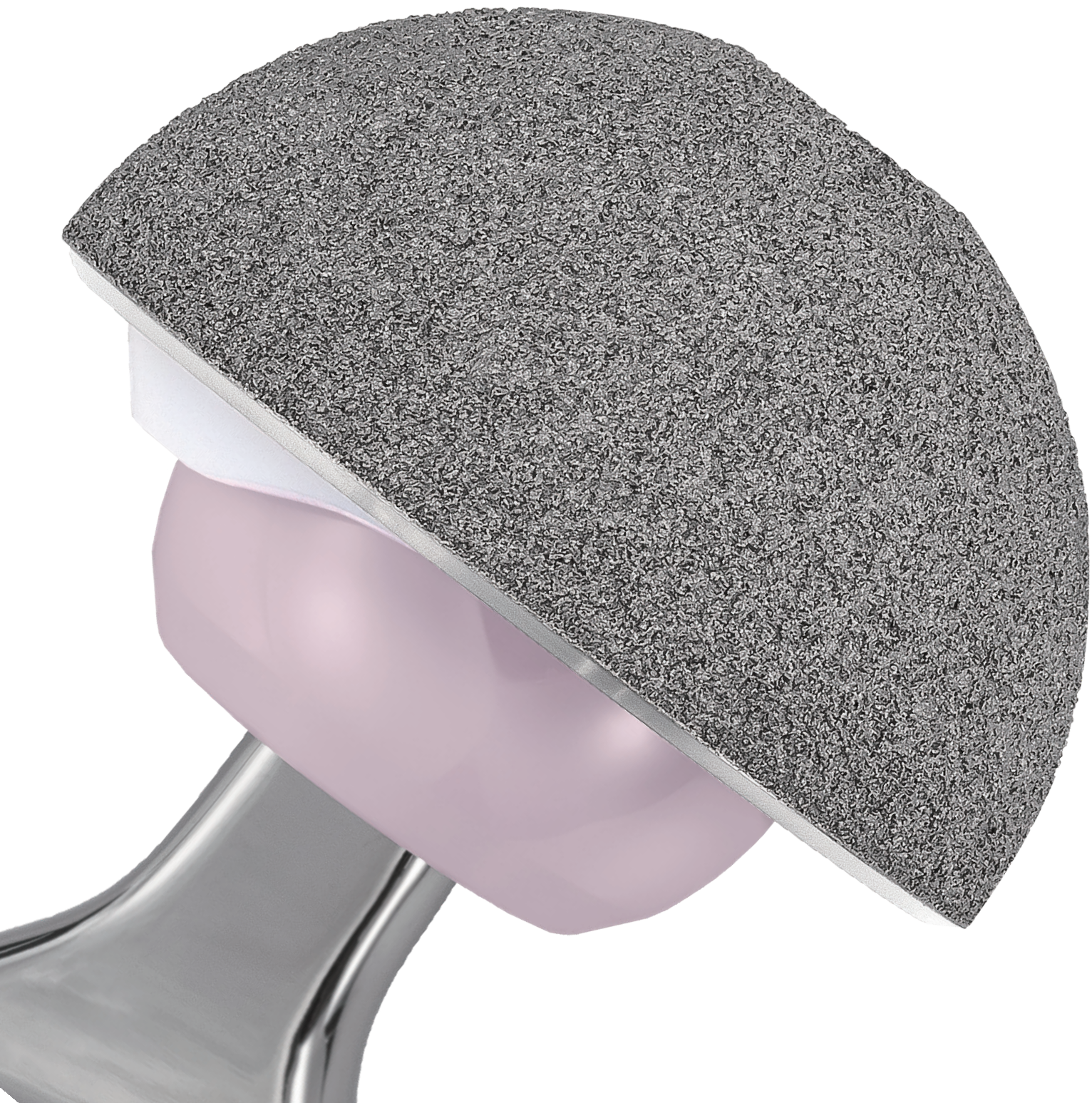

TPS Plus™ Surface Coating

U-Motion II⁺ Acetabular System utilizes the new Ti plasma spray plus (TPS Plus) coating technology to achieve a greater surface roughness (Ra).

The rough surface of implant provides scratch fit for initial stability and has been shown to reduce micromotion¹, both effects contribute to long term fixation of the implant^{1,2}.

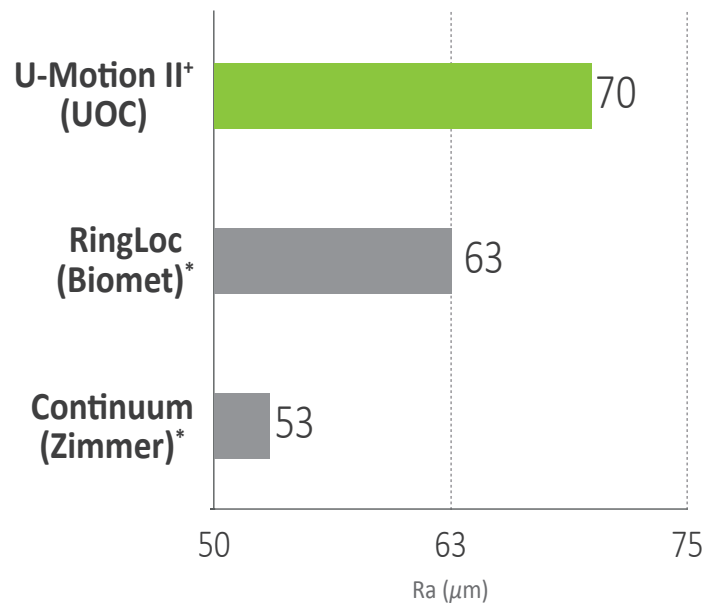

Increase Surface Roughness

**Data from RingLoc brochure. File at www.biomet.com*



Optimized Head-to-Cup Ratio Provides Excellent ROM

Larger head combined with 12/14 UTF taper stem system, allowing satisfied flexion range and greatly reduces impingement.




Wider Range of Motion




Simple Setting of Liner-Cup Assembly



Easy snap-in insertion

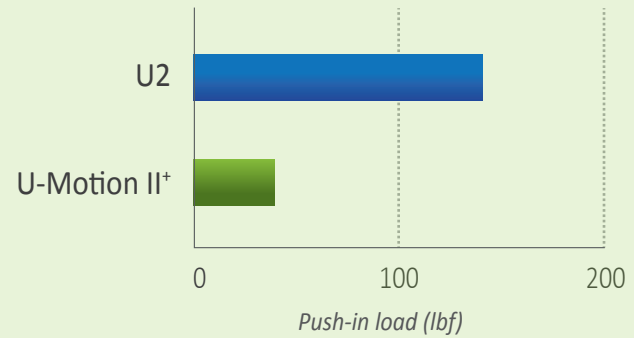

Enhanced locking strength

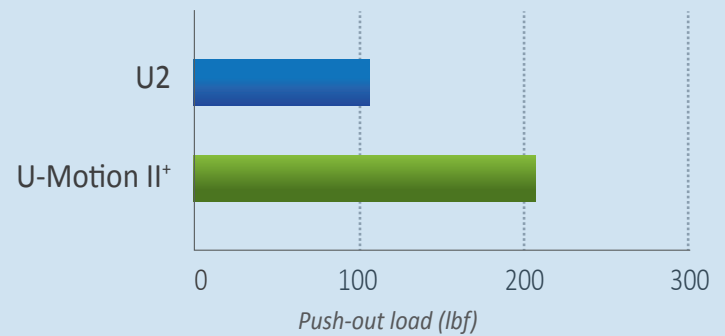

Data from U-Motion II+ testing report, UOC-RD-TR-12032



High Screw Angulation

Improving surgical flexibility without compromising cup fixation³.

Screw angulation




No-hole



Cluster-hole
with cover



Excellent Wear Performance
without the Risk of Oxidation

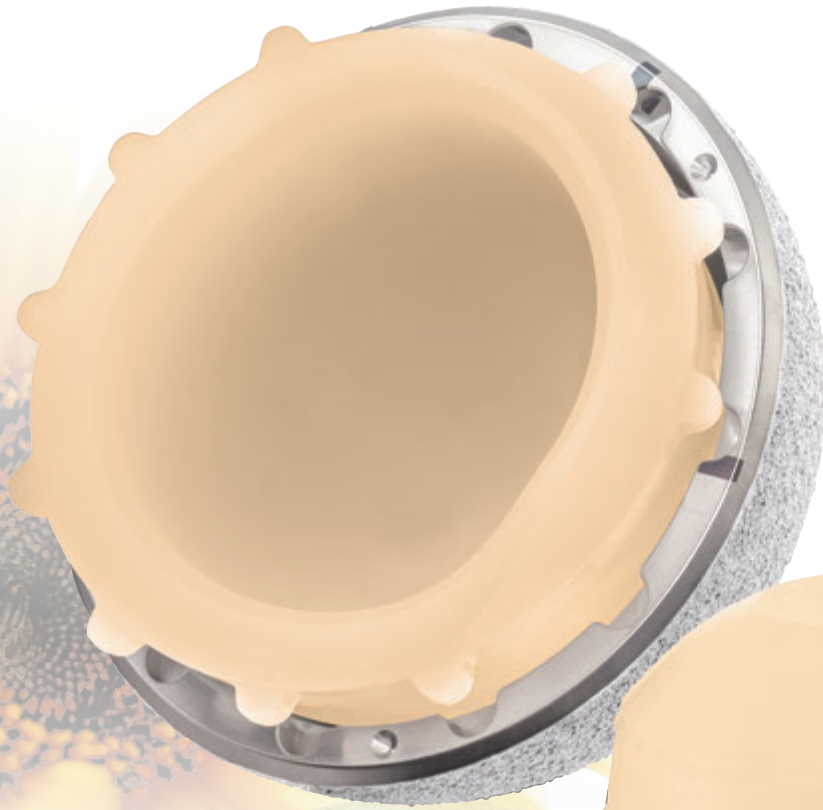


*E-XPE*TM

Vitamin E Stabilized Highly Cross-linked UHMWPE


Superior Oxidative Stability


Antioxidant Vitamin E provides long-term oxidative protection by stabilizing the free radicals⁴. Surface oxidative index of E-XPE showed no oxidation after in vitro accelerated aging test whereas the conventional highly cross-linked polyethylenes (HXLPE) do not eliminate oxidative potential⁵.

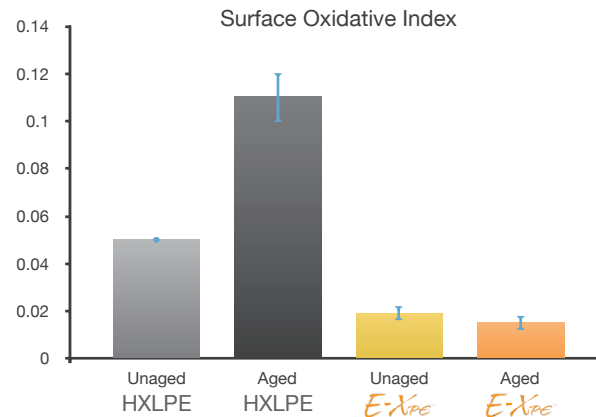

Extraordinary Wear Performance


125kGy radiation dosage optimizes the cross-link density and enhances the wear resistance of E-XPE over HXLPE products. E-XPE liner performed 75% reduction in gravimetric wear after artificial aging simulation corresponding to 10 years real-time usage^{6,7}.

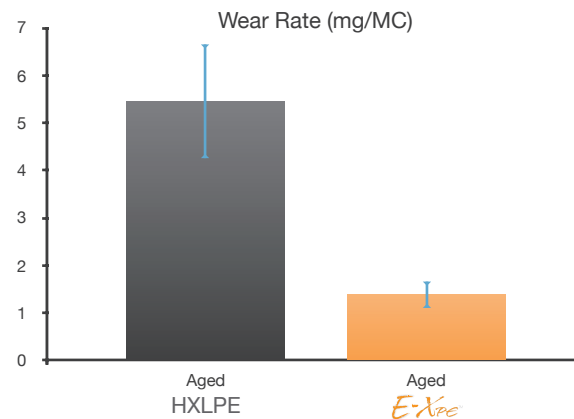

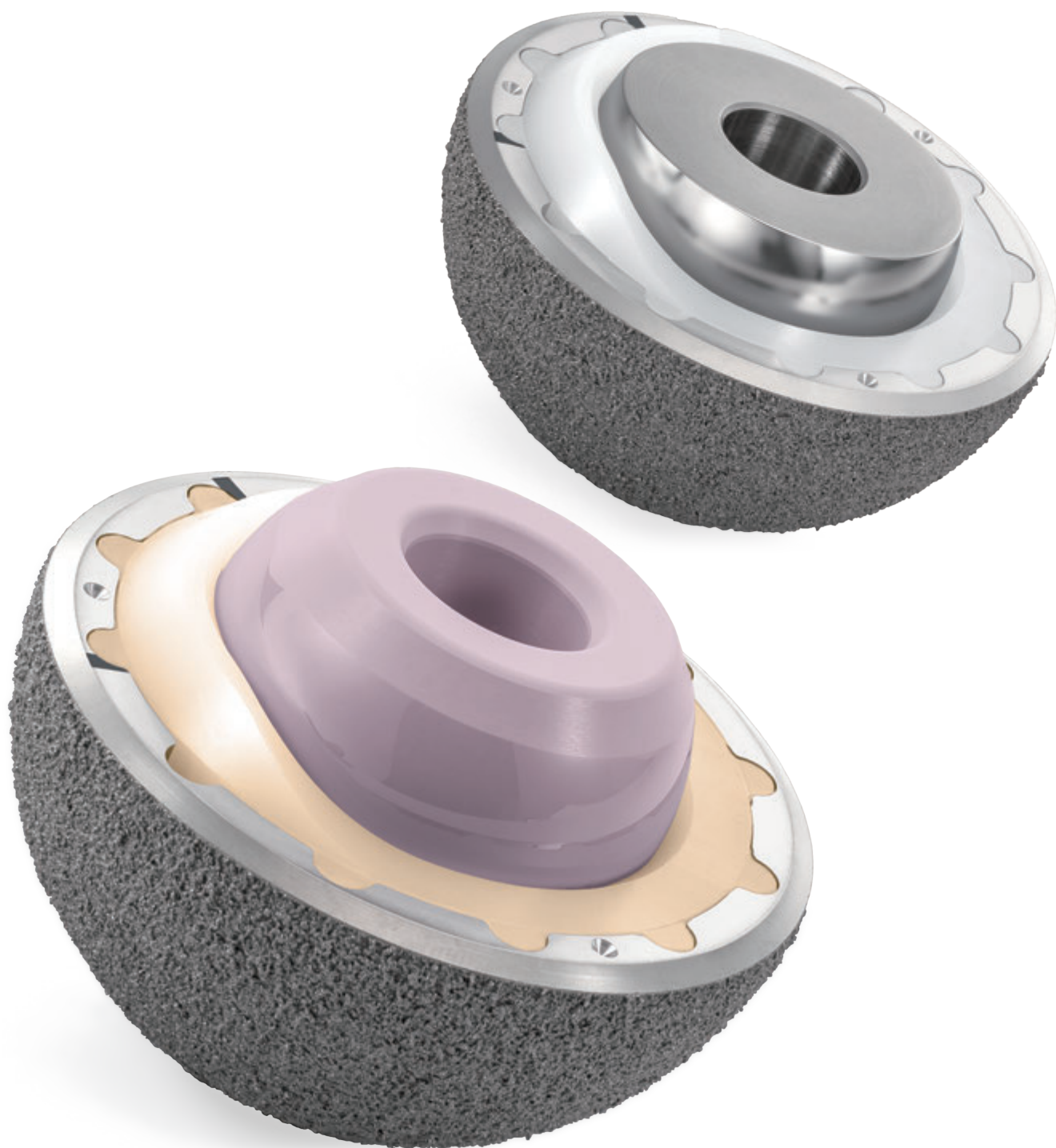


Bearing Options: Metal and Ceramic on E-XPE Liner

U-Motion II⁺ Acetabular System offers surgeons a wide range of bearing options, including both highly cross-linked polyethylene (XPE) and Vitamin E stabilized highly cross-linked polyethylene (E-XPE) liner that can accommodate both metal and ceramic heads.

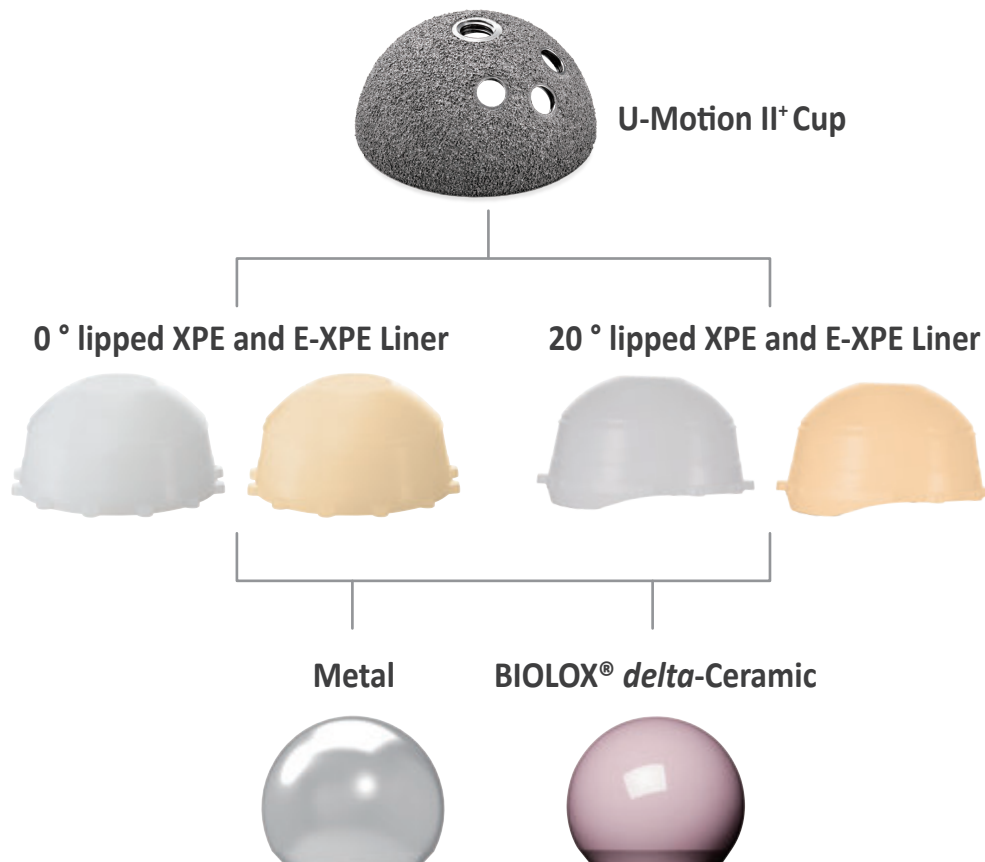

Order Information

★ *Optional Item*

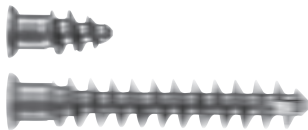

Ti Cancellous Screw, Ø 6.5 mm	
Length (mm)	Cat. No.
15	5206-1015
20	5206-1020
25	5206-1025
30	5206-1030
35	5206-1035
40	5206-1040
★ 45	5206-1045
★ 50	5206-1050



Screw Hole Covers	
Cat. No.	1306-1001

Screw hole covers are pre-assembled with multi-hole cup.



	Cluster-Hole Cup	No-Hole Cup
Cup Size (mm)	TPS Plus	TPS Plus
Ø 44	1306-3144	1306-3344
Ø 46	1306-3146	1306-3346
Ø 48	1306-3148	1306-3348
Ø 50	1306-3150	1306-3350
Ø 52	1306-3152	1306-3352
Ø 54	1306-3154	1306-3354
Ø 56	1306-3156	1306-3356
Ø 58	1306-3158	1306-3358
Ø 60	1306-3160	1306-3360
Ø 62	1306-3162	1306-3362
★ Ø 64	1306-3164	1306-3364
★ Ø 66	1306-3166	1306-3366
★ Ø 68	1306-3168	1306-3368
★ Ø 70	1306-3170	1306-3370

★ *Optional Item*



	XPE Liner		E-XPE Liner		Head Size	
Cup Size (mm)	Standard	20° lipped	Standard	20° lipped	CoCrMo	<i>delta</i> Ceramic
Ø 44	1406-3844	1406-5844	1406-7144	1406-7544	28 mm	28 mm
Ø 46	1406-3846	1406-5846	1406-7146	1406-7546		
Ø 48	1406-3248	1406-5248	1406-7248	1406-7648	32 mm	32 mm
Ø 50	1406-3250	1406-5250	1406-7250	1406-7650		
Ø 52	1406-3652	1406-5652	1406-7352	1406-7752	36 mm	36 mm
Ø 54	1406-3654	1406-5654	1406-7354	1406-7754		
Ø 56	1406-3656	1406-5656	1406-7356	1406-7756		40 mm
Ø 58	1406-3658	1406-5658	1406-7358	1406-7758		
Ø 60	1406-3660	1406-5660	1406-7360	1406-7760		
Ø 62	1406-3662	1406-5662	1406-7362	1406-7762		
★ Ø 64	1406-3664	1406-5664	1406-7364	1406-7764		
★ Ø 66	1406-3666	1406-5666	1406-7366	1406-7766		
★ Ø 68	1406-3668	1406-5668	1406-7368	1406-7768		
★ Ø 70	1406-3670	1406-5670	1406-7370	1406-7770		

★ *Optional Item*



Head						
Cup Size (mm)	CoCrMo			<i>delta</i> Ceramic		
Ø 44 Ø 46	1201-1028	-3	28 mm	1203-5028	-2.5	28 mm
	1201-1128	+ 0	28 mm	1203-5228	+ 1	28 mm
	1201-1228	+ 2.5	28 mm	1203-5428	+ 4	28 mm
	1201-1428	+ 5	28 mm			
	1201-1628	+ 7.5	28 mm			
	1201-1828	+ 10	28 mm			
Ø 48 Ø 50	1201-1032	-3	32 mm	1203-5032	-3	32 mm
	1201-1132	+ 0	32 mm	1203-5232	+ 1	32 mm
	1201-1232	+ 2.5	32 mm	1203-5432	+ 5	32 mm
	1201-1432	+ 5	32 mm	1203-5632	+ 8	32 mm
	1201-1632	+ 7.5	32 mm			
	1201-1832	+ 10	32 mm			
Ø 52 Ø 54				1203-5036	-3	36 mm
				1203-5236	+ 1	36 mm
				1203-5436	+ 5	36 mm
	1201-1036	-3	36 mm	1203-5636	+ 9	36 mm
Ø 56	1201-1136	+ 0	36 mm			
Ø 58	1201-1236	+ 2.5	36 mm			
Ø 60	1201-1436	+ 5	36 mm	1203-5040	-3	40 mm
Ø 62	1201-1636	+ 7.5	36 mm	1203-5240	+ 1	40 mm
★ Ø 64	1201-1836	+ 10	36 mm	1203-5440	+ 5	40 mm
★ Ø 66				1203-5640	+ 9	40 mm
★ Ø 68						
★ Ø 70						

*BIOLOX® *delta* is the registry trademark of Ceramtec AG.



Each Step
We Care



Each Step
We Care



References

1. Jasty M et al. In vivo skeletal responses to porous-surfaced implants subjected to small induced motions. J Bone Joint Surg Am. 1997 May;79(5):707-14.
2. Feighan JE et al. The influence of surface-blasting on the incorporation of titanium-alloy implants in a rabbit intramedullary model. J Bone Joint Surg Am. 1995 Sep;77(9):1380-95.
3. Hsu JT et al. Effects of screw eccentricity on the initial stability of the acetabular cup. Int Orthop. 2007; 31: 451-55.
4. Oral et al., The effect of alpha-tocopherol on the oxidation and free radical decay in irradiated UHMWPE. Biomaterials, 27(32): 5580–5587, 2006.
5. Data held on file. United Orthopedic Corporation.
6. Data held on file. United Orthopedic Corporation.
7. Schaffner and Zurbrugg, The effect of accelerated aging mechanical properties and oxidation of UHMWPE liners gamma-sterilized in nitrogen. 46th Annual Meeting, Orthopaedic Research Society, Orlando, pp. 545-562, 2000.



Each Step **We Care**

UOC USA Inc.

20 Fairbanks, Suite 173
Irvine, CA 92618
Tel: (949) 328-3366
Fax: (949) 328-3367

For our valued distribution partners, you can reach our customer service associates at customerservice@uocusa.com

For assistance on general and/or product related inquiries, please email us at sales@uocusa.com

Please refer to the product-specific package inserts for important information, including indications, contraindications, warnings, precautions, and potential adverse effects.



www.uocusa.com

Copyright © 2017 United Orthopedic Corp. All rights reserved.

0199-B-09-03-US